

MIDMARK®

Installation/Operation Instructions

9A179001/003 FIXED ARMBOARD

KIT COMPONENTS

Part Number	Description	Qty.
9A179001/002	Procedures Armboard	1
003-0654-00	Installation / Operation Instructions	1

REQUIRED TOOLS

None

INTRODUCTION


To hold an arm in a secure position while the table is in a flat position. Armboard has limited up, down and side movement.



WARNING


Mount 9A179001/002 Armboard on sidrail bar of exam table. Use caution when positioning Armboard and/or table. Do not use with any other accessories on same side of table. Always support Armboard assembly when changing positions on it. Assure all positions are solidly locked to prevent sudden movement during procedure.

INSTALLATION

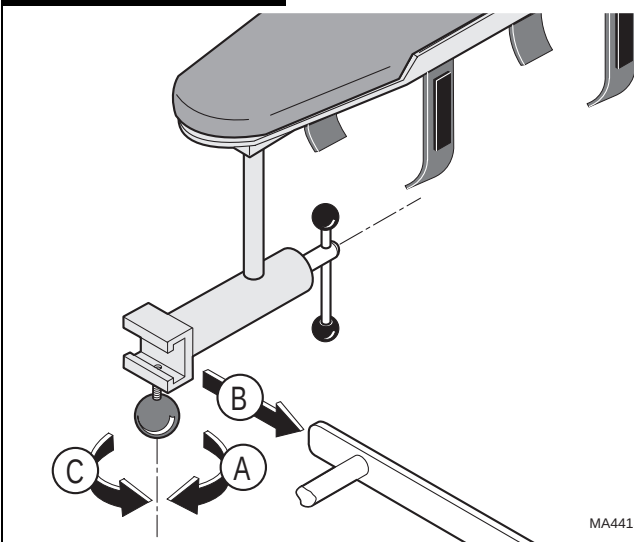

- A. Loosen knob.
- B. Slide swivel block onto side rail to desired position.
- C. Tighten knob.

OPERATION (HEIGHT & ROTATION)

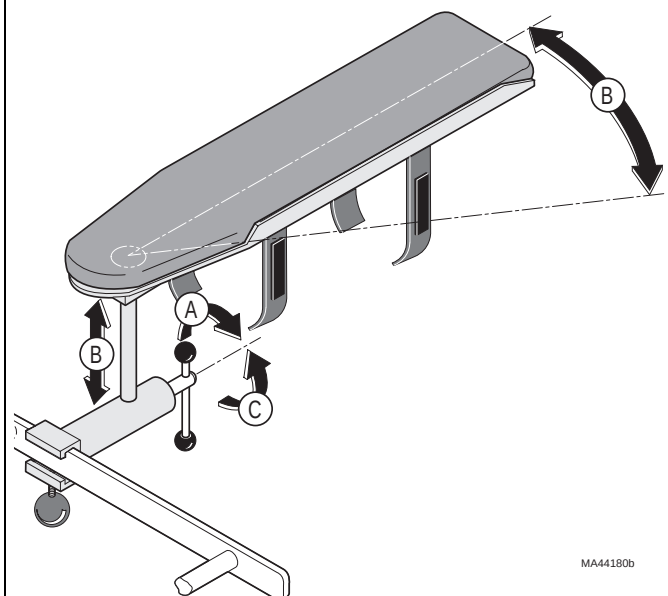

- A. Supporting Armboard, loosen T-handle.
- B. Raise or lower and rotate the Armboard to desired position.
- C. Tighten T-handle.

MIDMARK®

Midmark Corporation, Versailles, Ohio 45380 U.S.A.
(937) 526-3662 FAX (937) 526-5542

003-0654-00 Rev. B (4/98)